

## Technical information Serenys

63 69 75	39 45 51	22-29	52-45 (front) 38-53 (back)	102	117	40-44	50-60	-2°/+15°	-5°/+50°	19°/113°	120

Deviations in sizes and weight (1,5 cm and 1,5 kg) are possible

## General options



- Collapsible table (B12)



- Anti-tipping  
- Spiderwheels



- Side support



- Type 30

## The wheelchair must be adapted to you, not vice versa

The wheelchair has been developed as a modern aid which can be individually adapted. Aesthetic appearance and practical use go hand in hand. The wheelchair manufacturer Vermeiren has played an important role in this development.

Our basic principle is: the wheelchair must be adapted to you, not vice versa. This is why our wheelchairs have a modular design. All types of wheelchairs can be adapted, have numerous accessories and all options are mutually interchangeable. This is the result of our continuous research and more than 45 years of experience

## 3 Color combinations

Frame	Upholstery
Blue C20	Dartex Grey
Green C41	Dartex Green
Grey C29	Dartex Grey

## Your guarantee of quality

Regarding quality, the Serenys complies with the highest norms laid down by certified test institutes in most European countries.

## N.V. Vermeiren N.V.

Vermeirenplein 1-15  
B-2920 Kalmthout  
Belgium

Tel.: +32 (0) 3 620 20 20

Fax: +32 (0) 3 666 48 94


www.vermeiren.com • info@vermeiren.be



# VERMEIREN

Opening up new horizons in mobility

For more information, contact your local dealer

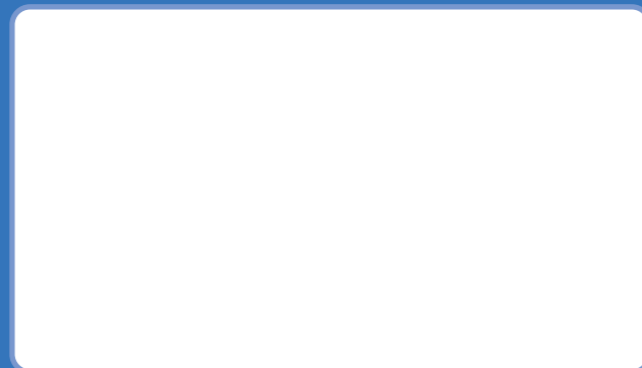

# Serenys

COMFORTABLE RELAXING

R.E.: N.V. Vermeiren N.V., Vermeirenplein 1/15 - 2920 Kalmthout - Belgium - 09/2005

Subject to technical alterations



# Serenys



Pushing handles adjustable in height.

Headrest adjustable in angle, height and depth.



The legrests are fully adjustable and the footrests with inclination will result in an ideal angle for the feet at any given time.



Removable back, seat and wheels ease the transport of the wheelchair.



Removable armrests which can be adjusted in height.



Back - high : 60 cm (to order)  
Back - low : 50 cm (standard version)



Standard removable armrests (not adjustable in height).

Seating height adjustable in 3 positions



Angle and height of the seating pillow is adjustable. Seating depth is possible in two sizes.



1. The legrests are adjustable in height.
2. The pillows of the legrests are collapsible to the back to ease transfers.
3. The pillows of the legrests are adjustable in height, width and depth.